

Best[®] nomos[®]
healthcare for everyone



Best[®] Sonalis[™]
ULTRASOUND IMAGING SYSTEM

A Breakthrough in Ultrasound Technology — Best[®] Sonalis[™] Imaging System

Featuring SimulView[™] Ultrasound Technology


The Best[®] Sonalis[™] Ultrasound Imaging System provides superior visualization of HDR, LDR, RF or Cryosurgical procedures. Our patented SimulView[™] Technology provides simultaneous "live" views of the prostate in both planes, thereby increasing treatment accuracy and precision.

Additional Features

- Patented probe design with simultaneous imaging of sagittal and transverse planes
- Longitudinal array provides for 140 mm length of view encompassing the bladder, prostate, and perineum
- Superior HD Image resolution for improved implant accuracy, speed and physician confidence level
- New system software re-designed with clinical need specifications in mind
- Advanced drawing and editing tools which include user-defined line widths and colors for fiducial and anatomical markers
- Advanced modular software design provides for future upgrade path via in-house and independently developed technologies
- User-friendly patient and image database uses familiar MS Windows folder and file workflow
- Multiple brachytherapy grids with user-definable grid symbols and colors
- Independent focal zones for sagittal and transverse images
- Independent frequency selection for sagittal and transverse Images
- On screen grid calibration

Advanced Keyboard Technology

Totally sealed, self-healing antibacterial keyboard with SensoFoil™ Technology

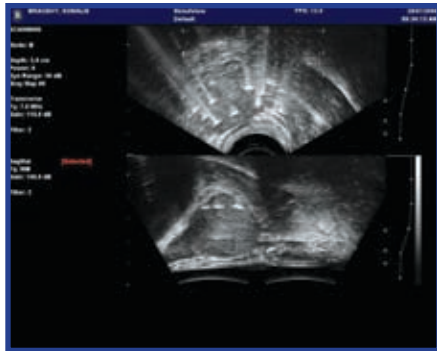

OUR SERVICE IS AS TRUSTWORTHY AS OUR PRODUCTS

We are committed to providing the highest level of customer training and support to ensure your practice runs smoothly without disruption to patient throughput.

SPECIALISTS ON CALL

Our specialists are on call Monday through Friday, 8 am to 5 pm EST. After-hours support is available via a paging system.

DEDICATED TO PROVIDING TREATMENT OPTIONS

At Best® NOMOS®, everything we do is designed to provide targeted, accurate and effective treatment for patients who are fighting cancer. By making our advanced technology available and affordable, we are helping more oncologists and clinicians treat more patients at more locations than ever before.

Best® NOMOS® is a member of TeamBest™
– a family of Best® Medical companies.

From brachytherapy seeds and equipment, dosimetry kits, phantoms, treatment planning systems, an array of medical/health physics equipment, and product remanufacturing/servicing, to radioactive sources, gamma teletherapy machines, cyclotron systems, and particle therapy treatment, TeamBest™ has it all!

Check out each company's site at www.teambest.com

Best® **nomos**®

1 Best Drive Pittsburgh, PA 15202 USA
Tel: 412 312 6700 800 70 NOMOS
Fax: 412 312 6701 www.nomos.com

Best® **medical international**

7643 Fullerton Road Springfield, VA 22153 USA
Tel: 703 451 2378 800 336 4970
Fax: 703 451 5228 www.teambest.com

BN/MB SONALIS – MSF-0001(R1) v.06.2016

© 2016 All Rights Reserved Best® NOMOS®. The specifications contained in this brochure are subject to change.

AFRICA | ASIA | EUROPE | LATIN AMERICA | MIDDLE EAST | NORTH AMERICA

